

Mr. Thirsty One-Step | High Volume Evacuator

Product Information Sheet



REF Mr. Thirsty® One Step | High Volume Evacuator

50Z988	Variety Pack 30 Count
50Z989LG	Large 50 Count
50Z989LG-B	Large 100 Count
50Z989MD	Medium 50 Count
50Z989MD-B	Medium 100 Count
50Z989SM	Small 50 Count
50Z989SM-B	Small 100 Count

Product Description

The Zirc Dental Products Mr. Thirsty One Step® is a high-volume evacuation medical device that connects to an HVE system during dental procedures. Mr. Thirsty One Step® features a polystyrene angled tube for connection to the HVE system and an adjustable thermoplastic tongue shield with integrated cheek flap and bite block.

Mr. Thirsty One Step® is a single-use, disposable device which is provided non-sterile. Mr. Thirsty One Step® devices are compatible with 7/16 in. diameter HVE valves.

Intended Use and Indication

Mr. Thirsty One Step® is intended for use with high volume evacuation for evacuation of fluids from the oral cavity while also retracting the tongue and cheek of the patient during dental procedures in which sterile devices are not required.

Specifications (size)

Mr. Thirsty One Step® is available in three sizes, each identified by a unique color (small- neon green, medium – neon purple and large – neon blue).

Warnings & Precautions

Caution: Federal (USA) law restricts this device to sale by or on order of a physician.



Do not use if packaging is damaged or unintentionally opened prior to use



Mr. Thirsty One Step® products are provided non-sterile



Do not reuse Mr. Thirsty One Step® High Volume Evacuators. Mr. Thirsty One Step® High Volume Evacuators are not designed to be cleaned, disinfected, or sterilized after use. Re-use could lead to microbial contamination and infection.




Keep product dry prior to use.



Keep product away from light sources prior to use.

* Users should cease use of Mr. Thirsty One Step® if it causes discomfort for the patient.

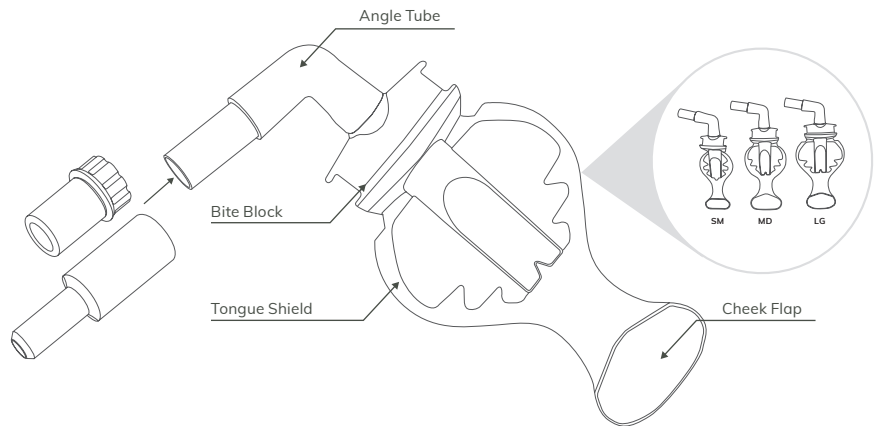

Use Specification

Intended Patient Population: General. No known contraindicated patient groups.

Part of the body interacted with: Oral cavity.

Intended User Profile: Dental professionals including Dentists, Dental Assistants, Dental Therapists, and Hygienists who have access to the current state-of-the-art of dentistry.

Intended Use Environment: State-of-the-art dentistry environment. No known restrictions on the use environment.

Medical Conditions Intended to Treat: There is no specific medical condition for which Mr. Thirsty® One-Step products are intended to treat. It is intended to be used, as needed, during dental procedures.

Operating Principle: Mr. Thirsty One Step® is designed to efficiently remove fluids from the oral cavity via direct connection to HVE systems. The integrated tongue-shield and bite block facilitate hands-free retraction and isolation.

Compatible Connections: Connect the Mr. Thirsty One-Step® angle tube directly to a 7/16-inch diameter HVE valve.

Use Instructions

1. Select a Mr. Thirsty One Step® High Volume Evacuator of the appropriate size. Mr. Thirsty One Step® is available in three sizes, each identified by a unique color:

- Small | Neon Green
- Medium | Neon Purple
- Large | Neon Blue

2. If adjustment is necessary, use scissors to trim around the perimeter of the tongue shield and cheek flap.

3. Fold cheek flap over Zirc logo.

4. Insert into mouth. Rotate to seat bite-block.

5. Connect to High Volume Evacuation Valve.

6. After use, dispose of Mr. Thirsty One Step® products as general waste. Mr. Thirsty One Step® High Volume Evacuators are **SINGLE USE**.



Zirc Dental Products, Inc.
3918 Hwy 55 SE
Buffalo, MN 55313
United States (US)
+1-763-682-6636



ifu.zirc.com

SYMBOLS GLOSSARY

 Indicates the product manufacturer.	 CAUTION: Federal law (U.S.) restricts this device to sale by or on the order of a physician.	 Product catalogue number.	 Product lot or batch code.	 Identifies the unique device identifier.
 Product date of manufacture.	 Keep dry.	 Keep away from direct sunlight.	 Do not use if package is damaged.	 Non-sterile.
 Do not re-use.	 Consult instructions for use or consult electronic instructions for use.			