
maxill

steri-sox™

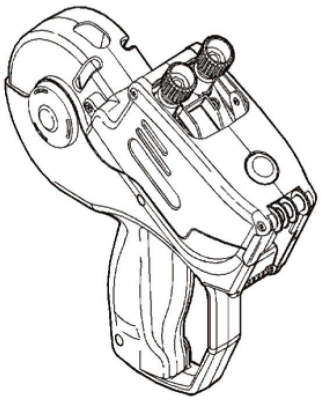


Instrument Reprocessing Labeler


Instructions For Use


Loading Labels

1. Push the release button and swing the cover open completely. You should hear a click. Don't force the cover past the click position.
2. For loading a new roll of supply, unwind the outer wrap of labels and liner.
3. Peel and discard the first 12 labels from the backing paper.
4. Place the label roll between the supply hubs.

   Do not thread label strip through any rollers.
5. Drape the supply across the labeler and center it over the peel bar. The supply should hang out 6-7 inches (15-16 cm) past the peel bar.
6. Close the cover.
7. Feed the label strip under the large roller and into the slot above the arrows. Hold the strip in place and squeeze the trigger until the label strip feeds through the rear exit.

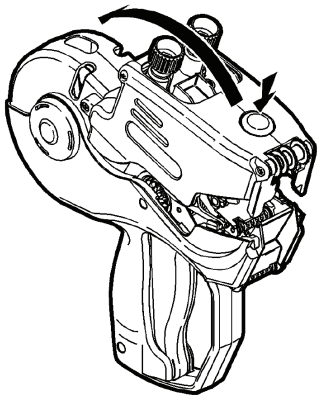


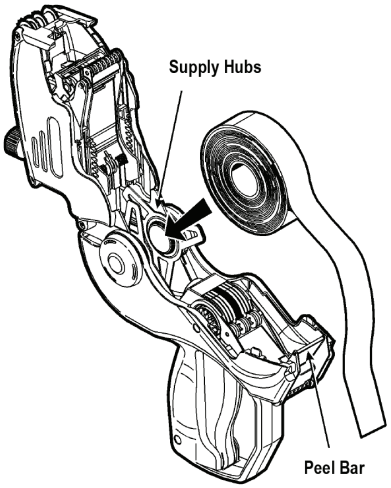

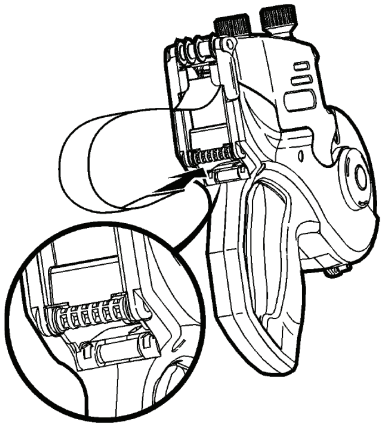

Setting Print Characters

1. Pull the knob out until the character pointers are on the desired position.
2. Turn the knob until the desired character with a white background appears between the pointers. To print a blank, select the blank white space. Characters with a black background should not appear in the window.
3. We do not recommend turning the knob past the white characters.
4. Push in the knob when you're finished.

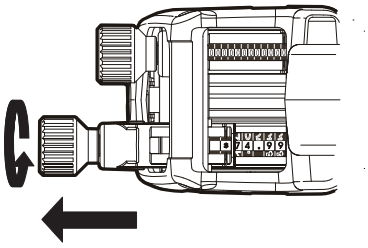

Using the Alpha Conversion chart

When using the middle alpha/numeric dial on the labeler what you read on the dial is not what is printed onto the label.

Use the Alpha Conversion chart located on the handle of the labeler to determine the letter you need to select to print the intended letter.

For example, if you read an *A* on the dial you will print a *T* on label. If you read a *T* you will print an *A*.

The letters on the middle row will read upside-down and backwards on the labeler but will print with the correct orientation.

Alpha Conversion


Read	or Print
A	T
B	U
C	V
D	W
E	X
F	Y
G	Z
H	1
J	2
K	3
L	4
M	5
N	6
P	7
Q	8
R	9
S	0

Special Instructions

Apply label to autoclave tape when using sterilization wrap.


Removing the Ink Roller

1. Open the cover completely.
2. Hold the ink roller on each end and push the release bar to release the ink roller.
3. Discard the used ink roller.

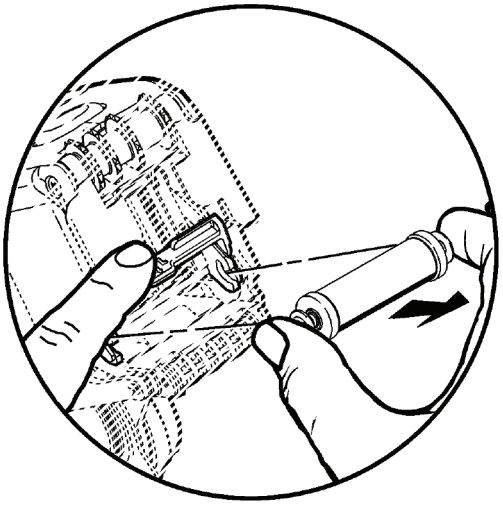


Installing the Ink Roller

1. Open the labeler cover completely.
2. Open the package and remove the ink roller.
3. Hold the new ink roller on each end.
4. Snap the ink roller into place.

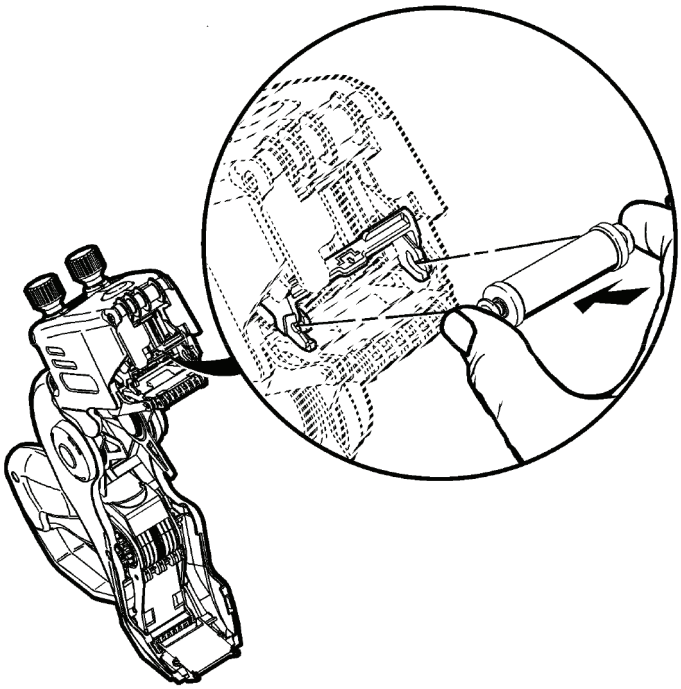

Clearing Label Jams, Cleaning, and Lubrication

1. Open the cover and remove the supply.
2. Remove any loose labels from the inside of the labeler.
3. Pull up on the deflector frame to expose the backing paper path.
4. Remove loose labels from under the deflector frame.
5. Snap the deflector frame back into place.
6. Use a commercial adhesive solvent to clean adhesive build-up
7. Reload and close the cover.

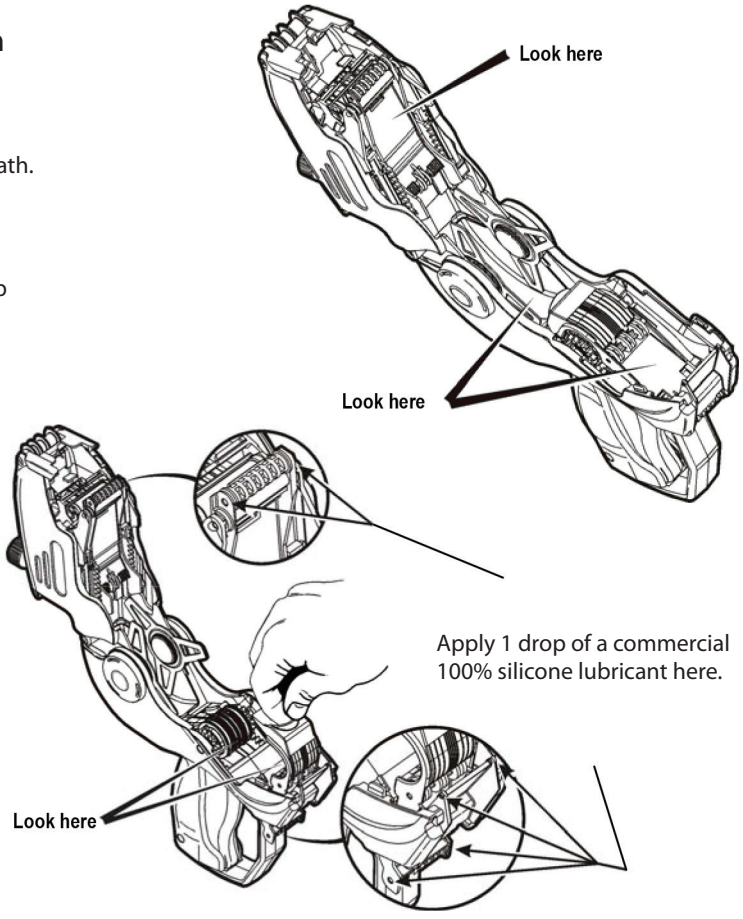

®maxill inc.
www.maxill.com 1-800-268-8633

USA: maxill inc.
Cortland, OH
44410 USA

Canada: maxill inc.
St. Thomas, ON
N5R 6C8 Canada

MP-REF-91140D
Revision 1: 2018-03-27