# Instructions For Use

maxill maxsėt




MIXED BERRY FLAVOUR

Item 53776

**TYPE I: Fast setting compliant with ISO 21563 specification.**

maxill maxsėt does not require special treatment before use.

**PURPOSE:** maxill maxsėt is a non-chromatic alginate with long dimensional stability specifically designed with the same optimal application specifications as maxill chrömamax. The orange colour allows an excellent contrast and tool to evaluate the quality of the impression.

maxill maxsėt alginate can be used for the preparation of:

- Preliminary models
- Study models
- Mobile prosthesis
- Removable and fixed prosthesis antagonist

**INSTRUCTIONS:**
Follow the powder/water ratio. Pour the powder into the water, and mix vigorously. Work in compliance with the declared times and chromatic phases. Measurements were taken at 23°C with deionized water. When using specific mechanical mixers for alginate-based impression materials, follow the manufacturer's instructions. Variations in powder/water ratio, temperature and/or water hardness, and speed and/or friction during mixing can cause changes in the working time and the initial setting time.

**WARNINGS:** For use, wear face mask, protective gloves, eye protection, and protective clothing. The powder is "Dust Free", however, cation is recommended to avoid inhalation. In case of contact with eyes wash thoroughly with water. After using the powder, close the bag immediately. Store in a cool and dry place. Use within the expiry date. Do not use on patients allergic to berry flavours.

Spaulding's Classifications

maxill inc.

**Canadian Head Office**
80 Elm Street
St. Thomas, ON N5R 6C8
Toll Free: 1-800-268-8633
Phone: 1-519-631-7370
Fax: 519-631-3388

**USA Head Office**
500 West Main Street
Cortland, OH 44410
Toll Free: 1-855-462-9455
Phone: 1-330-637-1176
Fax: 1-330-637-1154

Rev.0: 2025-01-09
MP-REF-53776D

maxill

# Instructions For Use


maxill maxsèt

Item 53776

**Made in Italy**



**20 mL water** (≈1 scoop) + **9 g maxsèt** (≈1 scoop)

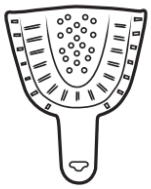

**40 mL water** (≈2 scoops) + **18 g maxsèt** (≈2 scoops)

**1**



**2**



**3**



**4**



**5**



MIXING TIME* **45s** — **40s**

WORKING TIME* **1m 25s**

INITIAL SETTING TIME* **1m 35s**

**6**





***ONLY* 30 s** *in the mouth**

(*) Measurements taken at 23°C with deionised water.

Spaulding's Classifications


**maxill inc.**

**Canadian Head Office**
80 Elm Street
St. Thomas, ON N5R 6C8
Toll Free: 1-800-268-8633
Phone: 1-519-631-7370
Fax: 519-631-3388

**USA Head Office**
500 West Main Street
Cortland, OH 44410
Toll Free: 1-855-462-9455
Phone: 1-330-637-1176
Fax: 1-330-637-1154



Rev.0: 2025-01-09

MP-REF-53776D