E-Z Jett® Cassettes

Product Demo Videos available online:
https://www.youtube.com/user/ZircDental

Zirc Company
3918 Hwy 55 S.E. | Buffalo, MN 55313 | USA
Toll Free: 800.328.3899 | Local: 763.682.6636
www.zirc.com


Zirc
Because Time Is Everything

Features:

- Made of durable, autoclavable, resin material
- Contains antimicrobial protection

Do not allow cassette to touch side walls or top/bottom of sterilizer

Squeeze tabs and lift to release cover

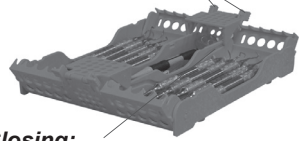Closing:
Pegs on cover should be in grooves on cassette


10-Place
(50Z930)

Squeeze tabs and lift to release cover

8-Place
(50Z927)

Squeeze tabs and lift to release cover

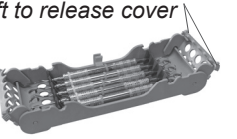5-Place
(50Z925)

Pull up cover to release

8-Place Slim
(50Z911)

Zirc Company

Processing of Resterilizable Medical Devices

Manufacturer/ Hersteller: Zirc Company, 3918 Hwy 55 SE, Buffalo, MN 55313 USA

1. Restriction of Reprocessing: Frequent reprocessing has marginal consequences for Zirc products. Normally, the end of durability of the product will be determined by visual degradation.
2. Cleaning: An Ultrasonic, Instrument Washer or Thermal Disinfector may be used. Follow manufacturer's directions for proper cleaning and use a mild alkaline detergent only. Rinse Zirc products with demineralized water to remove detergent or particles.
3. Packaging: Bag or wrap product using sterilization pouch or sterilization wrap.
4. Sterilization: Steam sterilize for at least 5 minutes at 273˚ F/ 134˚ C or 20 minutes at 250˚ F/ 121˚ C. Do not heat above 277˚ F/ 136˚ C. DO NOT dry heat sterilize.

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EC REP Loser & Co. GmbH Benzstr. 1-3, D51381 Leverkusen Germany +49-2171-706670